


OPRM1 Polyclonal Antibody

Product Details	
Size	100 µL
Species Reactivity	Human, Mouse, Rabbit, Rat
Published Species	Rat, Human, Mouse
Host/Isotype	Rabbit / IgG
Class	Polyclonal
Type	Antibody
Conjugate	Unconjugated
Immunogen	The antiserum was produced against a chemically synthesized peptide derived from an internal region of the human m-opioid receptor.
Form	Liquid
Purification	Antigen affinity chromatography
Storage buffer	Dulbecco's PBS, pH 7.3, with 1mg/mL BSA
Contains	0.05% sodium azide
Storage conditions	-20°C
RRID	AB_2533629

Applications	Tested Dilution	Publications
Western Blot (WB)	1:1,000	3 Publications
Immunohistochemistry (IHC)	-	1 Publication
Immunocytochemistry (ICC/IF)	1:250	1 Publication
Miscellaneous PubMed (Misc)	-	1 Publication

Product Images For OPRM1 Polyclonal Antibody

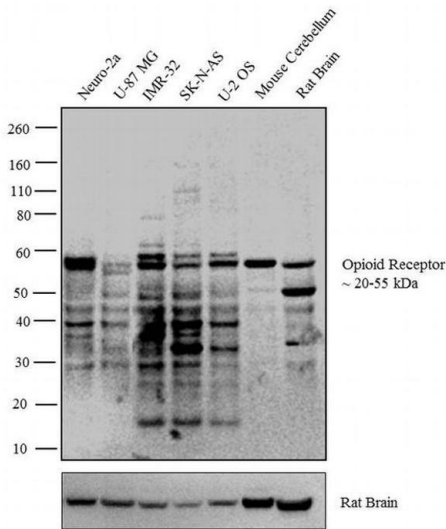

OPRM1 Antibody (44-308G) in WB



Western blot analysis was performed on membrane enriched extracts (30 µg lysate) of Neuro-2a (Lane 1), U-87 MG (Lane 2), IMR-32 (Lane 3), SK-N-AS (Lane 4), U-2 OS (Lane 5) and tissue extracts of Mouse Cerebellum (Lane 6) and Rat Brain (Lane 7). The blot was probed with Anti-Opioid Receptor Rabbit Polyclonal Antibody (Product # 44-308G, 1:250 dilution) and detected by chemiluminescence using Goat anti-Rabbit IgG (Heavy Chain) Superclonal™ Secondary Antibody, HRP conjugate (Product # A27036, 0.4 µg/mL, 1:2500 dilution). Opioid Receptor has multiple isoforms ranging from 20 to 55 kDa, which appears as a ladder form across the cell lines and tissues tested. Known quantity of protein samples were electrophoresed using Novex® NuPAGE and reg; 10% Bis-Tris gel (Product # NP0302BOX), XCell SureLock™ Electrophoresis System (Product # EI0002) and Novex® Sharp Pre-Stained Protein Standard (Product # LC5800). Resolved proteins were then transferred onto a nitrocellulose membrane with iBlot® 2 Dry Blotting System (Product # IB21001). The membrane was probed with the relevant primary and secondary Antibody following blocking with 5 % skimmed milk. Chemiluminescent detection was performed using Pierce™ ECL Western Blotting Substrate (Product # 32106).

OPRM1 Antibody (44-308G) in ICC/IF

Immunofluorescence analysis of OPIOID RECEPTOR was performed using 70% confluent log phase Neuro-2A cells. The cells were fixed with 4% paraformaldehyde for 10 minutes, permeabilized with 0.1% Triton™ X-100 for 10 minutes, and blocked with 1% BSA for 1 hour at room temperature. The cells were labeled with OPRM1 Rabbit Polyclonal Antibody (Product # 44-308G) at 1: 250 dilution in 0.1% BSA and incubated for 3 hours at room temperature and then labeled with Goat anti-Rabbit IgG (Heavy Chain) Superclonal™ Secondary Antibody, Alexa Fluor® 488 conjugate (Product # A27034) at a dilution of 1:2000 for 45 minutes at room temperature (Panel a: green). Nuclei (Panel b: blue) were stained with SlowFade® Gold Antifade Mountant with DAPI (Product # S36938). F-actin (Panel c: red) was stained with Rhodamine Phalloidin (Product # R415, 1: 300). Panel d represents the merged image showing membranous localization. Panel e shows the no primary antibody control. The images were captured at 60X magnification.

OPRM1 Antibody (44-308G) in WB

Western Blot. Extracts of rat brain lysates were resolved by SDS-PAGE on a 10% Tris-glycine gel and transferred to PVDF. The membrane was blocked with a 5% BSA-TBST buffer for one hour at room temperature and incubated with the μ-opioid receptor antibody for two hours at room temperature in a 1% BSA-TBST buffer. After washing, the membrane was incubated with goat F (ab')₂ anti-rabbit IgG HRP conjugate (Product # ALI4404) and signals were detected using the Pierce SuperSignal™ method

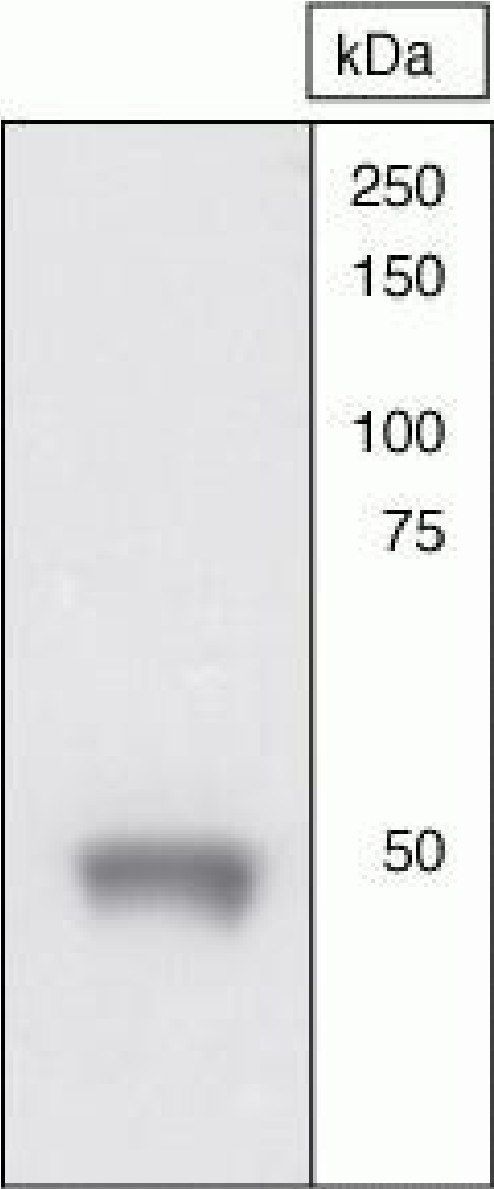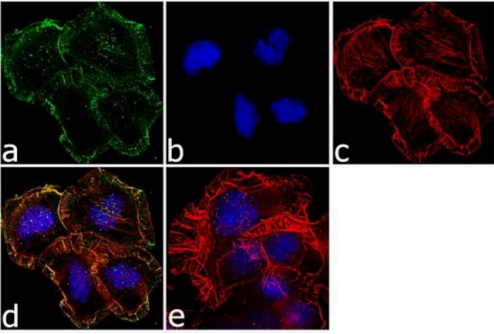

View more figures on thermofisher.cn

Western Blot (3)

<p>Neural regeneration research</p> <p>Effects of electroacupuncture on pain sensation in a rat model of hyperalgesia with nicotine dependence.</p> <p>"44-308G was used in Western Blotting to suggest that electroacupuncture treatment has positive analgesic effects on pain sensitivity caused by long-term chronic nicotine exposure."</p> <p>Authors: Wang SJ,Zhang YP,Candiotti KA</p>	<p>Year 2022</p> <p>Species Rat</p> <p>Dilution 1:300</p>
<p>Biochimica et biophysica acta</p> <p>Familial hemiplegic migraine type 1 mutations W1684R and V1696I alter G protein-mediated regulation of Ca(V)2.1 voltage-gated calcium channels.</p> <p>"44-308G was used in western blot to assess the effects of G protein-dependent modulation on mutations W684R and V696I which cause familial hemiplegic migraine type."</p> <p>Authors: Garza-López E,Sandoval A,González-Ramírez R,Gandini MA, Van den Maagdenberg A,De Waard M,Felix R</p>	<p>Year 2012</p> <p>Species Human</p> <p>Dilution 1:1000</p>

[View more WB references on thermofisher.cn](#)

Immunohistochemistry (1)

<p>The open orthopaedics journal</p> <p>Changes in midbrain pain receptor expression, gait and behavioral sensitivity in a rat model of radiculopathy.</p> <p>"44-308G was used in immunohistochemistry to assess gait and the expression of key pain receptors in the midbrain in a rodent model of radiculopathy."</p> <p>Authors: Hwang PY,Allen KD,Shamji MF,Jing L,Mata BA,Gabr MA,Huebner JL,Kraus VB,Richardson WJ,Setton LA</p>	<p>Year 2012</p> <p>Species Rat</p>
---	---

Immunocytochemistry (1)

<p>American journal of translational research</p> <p>Inhibition of mu-opioid receptor suppresses proliferation of hepatocellular carcinoma cells via CD147-p53-MAPK cascade signaling pathway.</p> <p>"44-308G was used in Immunocytochemistry to investigate the role of muopioid receptor (MOR) in the proliferation of Hepatocellular carcinoma (HCC) cell lines and the underlying mechanism."</p> <p>Authors: Zhang JJ,Song CG,Dai JM,Zhang XQ,Lin P,Li L,Yang XM,Chen ZN</p>	<p>Year 2021</p> <p>Species Human</p>
--	---

More applications with references on thermofisher.cn

Misc (1)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON-INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (i) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (ii) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (iii) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (iv) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.