# CD115 (c-fms) Monoclonal Antibody (AFS98), APC-eFluor™ 780, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Published Species	Fish, Mouse
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), APC-eFluor™ 780, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	AFS98
Conjugate	APC-eFluor™ 780
Excitation/Emission Max	756/785 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2573970

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.25 μg/test	25 Publications

# **Product Specific Information**

Description: The AFS98 monoclonal antibody reacts with the mouse CD115 molecule, a receptor for macrophage colony stimulating factor (M-CSF) or colony stimulating factor-1 (CSF-1). CD115 is expressed by monocyte, macrophage, osteoclast, and some epithelial cells. It is a 150 kDa c-fms gene product and belongs to immunoglobulin family. CSF-1 signaling through CSF-1R regulates the proliferation and differentiation of cells in the monocytic lineage.

Applications Reported: This AFS98 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This AFS98 antibody has been tested by flow cytometric analysis of mouse bone marrow-derived macrophages. This can be used at less than or equal to 0.25  $\mu$ g per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

APC-eFluor® 780 emits at 780 nm and is excited with the Red laser (633 nm). Please make sure that your instrument is capable of detecting this fluorochome.

Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.


Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222) (100 µL of cell sample + 100 µL of IC Fixation

Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 633-647 nm; Emission: 780 nm; Laser: Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.

# Product Images For CD115 (c-fms) Monoclonal Antibody (AFS98), APC-eFluor™ 780, eBioscience™

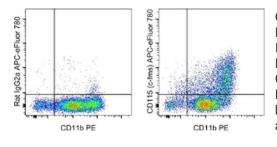

#### CD115 (c-fms) Antibody (47-1152-82) in Flow

Mouse bone marrow cells were cultured for 3 days with Mouse M-CSF Recombinant Protein (Product # 14-8983-80). Cells were stained with Anti-Mouse CD11b PE (Product # 12-0112-82) and 0.125 µg of Rat IgG2a K Isotype Control APC-eFluor® 780 (Product # 47-4321-82) (left) or 0.125 µg of Anti-Mouse CD115 (c-fms) APC-eFluor® 780 (right). Total viable cells, as determined by Fixable Viability Dye eFluor® 506 (Product # 65-0866-14), were used for analysis.

# View more figures on thermofisher.cn

□ 25 References	
Flow Cytometry (25)	
JCl insight Platelets exacerbate cardiovascular inflammation in a murine model of Kawasaki disease vasculitis.	<b>Year</b> 2023
"Published figure using CD115 (c-fms) monoclonal antibody (Product # 47-1152-82) in Flow Cytometry"	
Authors: Kocatürk B,Lee Y,Nosaka N,Abe M,Martinon D,Lane ME,Moreira D,Chen S,Fishbein MC,Porritt RA,Franklin BS,Noval Rivas M,Arditi M	
Frontiers in pharmacology	<b>Year</b> 2023
Pexidartinib synergize PD-1 antibody through inhibiting treg infiltration	
by reducing TAM-derived CCL22 in lung adenocarcinoma.	
"Published figure using CD115 (c-fms) monoclonal antibody (Product # 47-1152-82) in Flow Cytometry" Authors: Zhang W,Jiang X,Zou Y,Yuan L,Wang X	

View more Flow references on thermofisher.cn

# More applications with references on thermofisher.cn

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained Individuals. Unless otherwise stated in the Documentation, his warranty is package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained Individuals. Unless otherwise stated in the Documentation, this warranty are nor date of anyone of sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample furnished. Express OR INNER, DEPAR, REPLACE OR MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFERNING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAR, REPLACE COR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OPTION. THERE IS NO OBLICATION TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OPTION. THERE IS NO BULGATION TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagno

2