



CD41a Monoclonal Antibody (eBioMWReg30 (MWReg30)), Super Bright™ 436, eBioscience™

Product Details		
Size	100 μg	
Species Reactivity	Mouse	
Host/Isotype	Rat / IgG1, kappa	
Recommended Isotype Control	Rat IgG1 kappa Isotype Control (eBRG1), Super Bright™ 436, eBioscience™	
Class	Monoclonal	
Туре	Antibody	
Clone	eBioMWReg30 (MWReg30)	
Conjugate	Super Bright [™] 436	
Excitation/Emission Max	413/431 nm	
Form	Liquid	
Concentration	0.2 mg/mL	
Purification	Affinity chromatography	
Storage buffer	PBS, pH 7.2, with BSA	
Contains	0.09% sodium azide	
Storage conditions	4° C, store in dark, DO NOT FREEZE!	
RRID	AB_2784804	

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	1 Publication
Immunocytochemistry (ICC/IF)	-	1 Publication
Flow Cytometry (Flow)	1.0 µg/test	4 Publications

Product Specific Information

Description: The eBioMWReg30 monoclonal antibody reacts with mouse CD41 (fibrinogen receptor, gpllb, integrin alpha IIb). While initially thought to be expressed exclusively on the surface of platelets and megakaryocytes, it has been demonstrated that CD41 is also expressed on hematopoietic progenitors in the embryo, fetus and adult. CD41 associates with CD61 (gpllla, integrin beta III) to form a receptor which plays a major role in platelet function, including binding of several adhesion molecules such as fibrinogen, fibronectin and vitronectin.

Recently, the SLAM-family markers, CD48 and CD150 have been used to reliably identify hematopoietic stem cells (HSC). Specifically, it was found that CD150+CD48- bone marrow cells were highly efficient in their ability to confer long-term multi-lineage reconstitution in irradiated mice. Furthermore, the efficiency of reconstitution was enhanced when HSCs were further enriched through the exclusion of CD41+ cells. Thus, the use of CD150+CD48-CD41- as an expression profile efficiently identifies hematopoietic stem cells.

Applications Reported: This eBioMWReg30 (MWReg30) antibody has been reported for use in flow cytometric analysis.

Applications Tested: This eBioMWReg30 antibody has been tested by flow cytometric analysis of mouse platelets. This may be used at less than or equal to 1.0 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is

recommended that the antibody be carefully titrated for optimal performance in the assay of interest.


Super Bright 436 can be excited with the violet laser line (405 nm) and emits at 436 nm. We recommend using a 450/50 bandpass filter, or equivalent. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.

Excitation: 405 nm; Emission: 436 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD41a Monoclonal Antibody (eBioMWReg30 (MWReg30)), Super Bright™ 436, eBioscience™

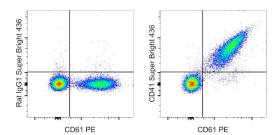

CD41a Antibody (62-0411-82) in Flow

C57BL/6 mouse platelets were stained with CD61 Monoclonal Antibody, PE (Product # 12-0611-82) and 0.5 µg of Rat IgG1 kappa Isotype Control, Super Bright 436 (Product # 62-4301-82) (left) or 0.5 µg of CD41 Monoclonal Antibody, Super Bright 436 (right).

View more figures on thermofisher.cn

□ 6 References

Immunohistochemistry (1)

British journal of haematology

All-trans-retinoic acid shifts Th1 towards Th2 cell differentiation by targeting NFAT1 signalling to ameliorate immune-mediated aplastic anaemia.

Year 2020

"Published figure using CD41a monoclonal antibody (Product # 62-0411-82) in Immunohistochemistry" Authors: Tang D,Liu S,Sun H,Qin X,Zhou N,Zheng W,Zhang M,Zhou H,Tuersunayi A,Duan C,Chen J

Immunocytochemistry (1)

Scientific reports

EphrinB2 regulates the emergence of a hemogenic endothelium from the aorta.

the aorta.
"Published figure using CD41a monoclonal antibody (Product # 62-0411-82) in Immunocytochemistry"

Year 2016

Flow Cytometry (4)

The Journal of biological chemistry

Authors: Chen II, Caprioli A, Ohnuki H, Kwak H, Porcher C, Tosato G

SNAP23 is essential for platelet and mast cell development and required in connective tissue mast cells for anaphylaxis.

"Published figure using CD41a monoclonal antibody (Product # 62-0411-82) in Flow Cytometry"

Authors: Cardenas RA,Gonzalez R,Sanchez E,Ramos MA,Cardenas EI,Rodarte AI,Alcazar-Felix RJ,Isaza A,Burns AR, Heidelberger R,Adachi R

Year 2021

Cell

Hematopoietic Stem Cells Count and Remember Self-Renewal Divisions.

"Published figure using CD41a monoclonal antibody (Product # 62-0411-82) in Flow Cytometry" Authors: Bernitz JM,Kim HS,MacArthur B,Sieburg H,Moore K

Year 2016

View more Flow references on thermofisher.cn

More applications with references on thermofisher.cn

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production of specifications and/or accompanying package inserts ("Documentation"), No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is more in a policy of the provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty does not extend to anyone other than the Buyer. Any mode of a sample in merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS OURING THE WARRANTY PERIOD IS LIMITED, CERPLIAR, PEPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS, AT SELLER'S SOLE OFFICIAL THE IS NO OBLIGATION TO REPLIAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS, AT SELLER'S SOLE OFFICIAL THE IS NO DELICATION TO REPLIAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS, AT SELLER'S SOLE OFFICIAL THE NON-CONFORMING PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vit