Performanc guarenteed

CD62L (L-Selectin) Monoclonal Antibody (MEL-14), Super Bright[™] 436, eBioscience[™]

Product Details

i louuci Detalis	
Size	100 µg
Species Reactivity	Mouse
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), Super Bright™ 436, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	MEL-14
Conjugate	Super Bright™ 436
Excitation/Emission Max	413/431 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2762739

Applications	Tested Dilution	Publications
Immunocytochemistry (ICC/IF)	-	1 Publication
Flow Cytometry (Flow)	0.25 µg/test	25 Publications

Product Specific Information

Description: The MEL-14 monoclonal antibody reacts with mouse CD62L, a 76 kDa member of the selectin family. CD62L is expressed by neutrophils, monocytes, and subsets of T, B, and NK cells and binds a number of glycosylated, fucosylated, sulfated sialylated glycoproteins including CD34, glycam-1 and MAdCam-1. These interactions mediate rolling of lymphocytes on activated endothelium at the sites of inflammation and homing of cells to the high endothelial venules (HEV) of peripheral lymphoid tissues.

Applications Reported: This MEL-14 antibody has been reported for use in flow cytometric analysis.


Applications Tested: This MEL-14 antibody has been tested by flow cytometric analysis of mouse splenocytes. This may be used at less than or equal to 0.25 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Super Bright 436 can be excited with the violet laser line (405 nm) and emits at 436 nm. We recommend using a 450/50 bandpass filter, or equivalent. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information. Excitation: 405 nm; Emission: 436 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD62L (L-Selectin) Monoclonal Antibody (MEL-14), Super Bright[™] 436, eBioscience[™]

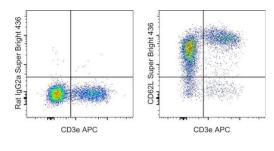

CD62L (L-Selectin) Antibody (62-0621-82) in Flow

C57BL/6 mouse splenocytes were stained with CD3e Monoclonal Antibody, APC (Product # 17-0031-82) and 0.125 μ g of Rat IgG2a kappa Isotype Control, Super Bright 436 (Product # 62-4321-82) (left) or 0.125 μ g of CD62L Monoclonal Antibody, Super Bright 436 (right). Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.cn

□ 26 References

Immunocytochemistry (1)

Cell communication and signaling : CCS Aiduqing formula inhibits breast cancer metastasis by suppressing TAM /CXCL1-induced Treg differentiation and infiltration.	Year 2021
"Published figure using CD62L (L-Selectin) monoclonal antibody (Product # 62-0621-82) in Immunocytochemistry"	
Authors: Li J,Wang S,Wang N,Zheng Y,Yang B,Wang X,Zhang J,Pan B,Wang Z	

Flow Cytometry (25)

NID Loss a store

Bioactive materials	Year
A platinum@polymer-catechol nanobraker enables radio-	2023
immunotherapy for crippling melanoma tumorigenesis, angiogenesis, and radioresistance.	
"Published figure using CD62L (L-Selectin) monoclonal antibody (Product # 62-0621-82) in Flow Cytometry"	
Authors: Li W,Yan J,Tian H,Li B,Wang G,Sang W,Zhang Z,Zhang X,Dai Y	

NPJ vaccines	2022	
Amphiphile-CpG vaccination induces potent lymph node activation and COVID-19 immunity in mice and non-human primates.		
"Published figure using CD62L (L-Selectin) monoclonal antibody (Product # 62-0621-82) in Flow Cytometry"		
Authors: Seenappa LM,Jakubowski A,Steinbuck MP,Palmer E,Haqq CM,Carter C,Fontenot J,Villinger F,McNeil LK, DeMuth PC		

View more Flow references on thermofisher.cn

More applications with references on thermofisher.cn

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/second by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Online work multiple and quality of goods and does not represent that any Product will conform to such model or sample furnitised to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample furnitised to Buyer is morely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample furnitised to Buyer is morely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample furnitised. EXERCENCE OR NON INFERINGENENT. BUYER'S EXCUSIVE ENERGIP FOR NON-CONFIGNING PRODUCTS DURING THE WARRANTES, EXPRESS OR INFENDIOS TO UNITED TEACHING TO REPLAR. REFLACE OR REFLIND FOR THE NON-CONFORNING PRODUCTS SOLE OF TOR NON-CONFIGNING PRODUCTS SOLE OF TOR UNATION TO FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER', (III) USE OF THE PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER', (III) USE OF THE PRODUCTS. UNITED STATE WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product is intended for research only and is not to be used for any ot

2