Performance guarenteed

NOX2 Polyclonal Antibody

Product Details

Size	100 µg
Species Reactivity	Human, Mouse, Rat
Published Species	Rat
Host/Isotype	Rabbit / IgG
Class	Polyclonal
Туре	Antibody
Conjugate	Unconjugated
Immunogen	Rabbit polyclonal NOX2 antibody was raised against a 15 amino acid peptide near the amino terminus of human NOX2. The immunogen is located within amino acids 140 - 190 of NOX2.
Form	Liquid
Concentration	1 mg/mL
Purification	Antigen affinity chromatography
Storage buffer	PBS
Contains	0.02% sodium azide
Storage conditions	Store at 4°C short term. For long term storage, store at -20°C, avoiding freeze/thaw cycles.
RRID	AB 2890479


Applications	Tested Dilution	Publications
Western Blot (WB)	1-2 µg/mL	-
Immunohistochemistry (IHC)	-	1 Publication
Immunohistochemistry (Paraffin) (IHC (P))	2 μg/mL	-

Product Specific Information

Positive controls that may be used include: Rat Brain Tissue Lysate

1

Product Images For NOX2 Polyclonal Antibody

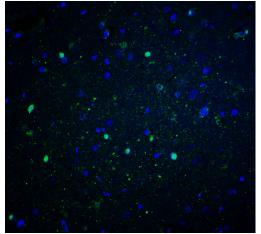


NOX2 Antibody (PA5-114489) in IHC (P)

Immunofluorescence of NOX2 in rat brain tissue with NOX2 Polyclonal Antibody (Product # PA5-114489) at 20 µg/mL. Green: NOX2 Blue: DAPI staining

NOX2 Antibody (PA5-114489) in IHC (P)

Immunohistochemistry of NOX2 in rat brain tissue cells with NOX2 Polyclonal Antibody (Product # PA5-114489) at 2 µg/mL.

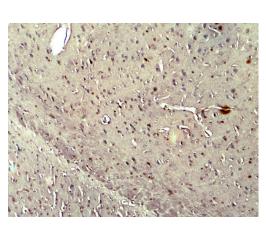

NOX2 Antibody (PA5-114489) in IHC (P) Immunohistochemistry of NOX2 in rat brain tissue with NOX2 Polyclonal Antibody (Product # PA5-114489) at 5 µg/mL.

View more figures on thermofisher.cn

□ 1 Reference

Immunohistochemistry (1)

Cells	Year
Hyperoxaluria Induces Endothelial Dysfunction in Preglomerular	2022
Arteries: Involvement of Oxidative Stress.	Species
"PA5-114489 was used in Immunohistochemistry-immunofluorescence to characterize urolithiasis-related endothelial dysfunction in a hyperoxaluria rat model of renal lithiasis."	Rat
Authors: Saenz-Medina J,Muñoz M,Rodriguez C,Contreras C,Sánchez A,Coronado MJ,Ramil E,Santos M,Carballido J, Prieto D	

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to oneyear from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLED WARRANTES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEEV FOR NON-CONFORMING PRODUCTS OURING THE WARRANTY FENDO IS LIMITED TO REPAR, REFLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS) AT SELLER'S SOLE OPTION. THERE JS NO BULGATION TO REPAR, REFLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS UNRING THE WARRANTY FENDO IS LIMITED TO REPAR, REFLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS) AT SELLER'S SOLE OPTION. THERE JS NO BULGATION TO REFAR, REFLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS UNRING THE WARRANT PENDO IS LIMITED TO REPAR, REFLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS) AT SELLER'S SOLE OPTION. THERE JS NO BULGATION TO REFAR, REFLACE OR REFUND FOR THE RODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJELIER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. UNRIGHT THE OPTION TO IN the documentation accompanying the Pro