Performanc guarenteed

CD8a Monoclonal Antibody (53-6.7), Biotin, eBioscience™

Product Details

| Size | 100 µg |
|--------------------------------|---|
| Species Reactivity | Mouse |
| Published Species | Fruit fly, Hamster, Mouse, Human |
| Host/Isotype | Rat / IgG2a, kappa |
| Recommended Isotype Control | Rat IgG2a kappa Isotype Control (eBR2a), Biotin, eBioscience™ |
| Class | Monoclonal |
| Туре | Antibody |
| Clone | 53-6.7 |
| Conjugate | Biotin |
| Form | Liquid |
| Concentration | 0.5 mg/mL |
| Purification | Affinity chromatography |
| Storage buffer | PBS, pH 7.2 |
| Contains | 0.09% sodium azide |
| Storage conditions | 4° C, store in dark, DO NOT FREEZE! |
| RRID | AB_466346 |
| | |

| Applications | Tested Dilution | Publications |
|---|-----------------|------------------|
| Western Blot (WB) | - | 2 Publications |
| Immunohistochemistry (IHC) | - | 12 Publications |
| Immunohistochemistry (Frozen) (IHC (F)) | Assay-Dependent | - |
| Immunocytochemistry (ICC/IF) | - | 6 Publications |
| Flow Cytometry (Flow) | 0.5 µg/test | 141 Publications |
| Immunoprecipitation (IP) | - | 1 Publication |
| Functional Assay (FN) | - | 3 Publications |
| Miscellaneous PubMed (Misc) | - | 8 Publications |

Product Specific Information

Description: The 53-6.7 monoclonal antibody reacts with the mouse CD8a molecule. CD8a is an approximately 32-34 kDa cell surface receptor expressed either as a heterodimer with the CD8 beta chain (CD8 alpha beta) or as a homodimer (CD8 alpha alpha). A majority of thymocytes and a subpopulation of mature alpha beta TCR T cells express CD8 alpha beta while gamma delta TCR T cells, a subpopulation of intestinal intraepithelial lymphocytes (IELs) and dendritic cells express CD8 alpha alpha. CD8 binds to MHC class I and through its association with protein tyrosine kinase p56lck plays a role in T cell development and activation of mature T cells.

Applications Reported: The 53-6.7 antibody has been reported for use in flow cytometric analysis and immunohistochemistry.


Applications Tested: The 53-6.7 antibody has been tested by flow cytometric analysis of mouse thymocytes and splenocytes. This can be used at less than or equal to 0.5 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It

1

is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For CD8a Monoclonal Antibody (53-6.7), Biotin, eBioscience™

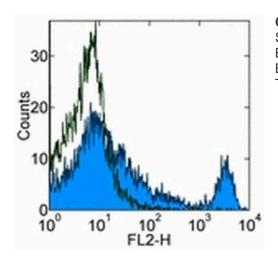

CD8a Antibody (13-0081-82) in Flow

Staining of C57BL/6 splenocytes with 0.25 µg of Rat IgG2a K Isotype Control Biotin (Product # 13-4321-82) (open histogram) or 0.25 µg of Anti-Mouse CD8a Biotin (filled histogram) followed by Streptavidin PE (Product # 12-4317-87). Total viable cells were used for analysis.

View more figures on thermofisher.cn

Western Blot (2)

| eLife | Year |
|---|------|
| Hemozoin produced by mammals confers heme tolerance. | 2019 |
| Published figure using CD8a monoclonal antibody (Product # 13-0081-82) in Western Blot" | |
| Authors: Pek RH,Yuan X,Rietzschel N,Zhang J,Jackson L,Nishibori E,Ribeiro A,Simmons W,Jagadeesh J,Sugimoto H, Alam MZ,Garrett L,Haldar M,Ralle M,Phillips JD,Bodine DM,Hamza I | |
| Cardiovascular research | Year |
| Genetic deletion of chemokine receptor Ccr7 exacerbates atherogenesis in ApoE-deficient mice. | 2013 |
| 'Published figure using CD8a monoclonal antibody (Product # 13-0081-82) in Flow Cytometry" Authors: Wan W,Lionakis MS,Liu Q,Roffê E,Murphy PM | |
| nmunohistochemistry (12) | |
| Scientific reports | Year |
| Therapeutic potential of highly functional codon-optimized | 2022 |
| microutrophin for muscle-specific expression. | |
| Published figure using CD8a monoclonal antibody (Product # 13-0081-82) in Immunohistochemistry" | |
| Authors: Starikova AV,Skopenkova VV,Polikarpova AV,Reshetov DA,Vassilieva SG,Velyaev OA,Shmidt AA,Savchenko IM,Soldatov VO,Egorova TV,Bardina MV | |
| Cancers | Year |
| TIAM2S Mediates Serotonin Homeostasis and Provokes a Pro- | 2020 |
| Inflammatory Immune Microenvironment Permissive for Colorectal | |
| Tumorigenesis. | |
| Published figure using CD8a monoclonal antibody (Product # 13-0081-82) in Immunohistochemistry" | |
| | |

More applications with references on thermofisher.cn

ICC/IF (6) Flow (141) IP (1) FN (3) Misc (8)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. No OTHER MARAINTES, EXPRESS or RMPLED, ARRAINTES, EXPRESS FOR ANY PARTICULAR PURPRES, GRANTE DINCES, CON NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS OURING THE WARRANTY PERIOD IS UNITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS) AS THER SOL DOTION. THERE IS NO OBLICATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS). USE OFTION. THERE IS NO OBLICATION TO REPAIR, REPLACE OR REFUND FOR THE RODUCTS. UNITED TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS). A MANNER FOR WHICH THEY WERE NOT DEFORMENT, BUYER, (III) MISUE, FUND FOR THE NON-CONFORMING PRODUCTS. UNITED TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS. UNITED TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS. UNITED TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS. UNITED TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS. UNITED TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS. UNITED TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS. UNITED TO REPAIR, REPLACE OR REFUND FOR THE NON-C