



CD366 (TIM3) Monoclonal Antibody (F38-2E2), PE-Cyanine5, eBioscience™

100 Tests
Human
Mouse / IgG1, kappa
Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), PE-Cyanine5, eBioscience™
Monoclonal
Antibody
F38-2E2
PE-Cyanine5
568/666 nm
Liquid
5 μL/Test
Affinity chromatography
PBS, pH 7.2, with 0.2% BSA
0.09% sodium azide
4° C, store in dark, DO NOT FREEZE!
AB_2802205
F F 0 4

Applications	Tested Dilution	Publications
Western Blot (WB)	-	1 Publication
Flow Cytometry (Flow)	5 μL (0.125 μg)/test	4 Publications

Product Specific Information

Description: This F38-2E2 monoclonal antibody reacts with human CD366, also known as T cell immunoglobulin and mucin domain-containing protein 3 (TIM3) or HAVCR2. This cell surface receptor is expressed on activated CD4+ T cell subsets (e.g. Th1, Th17, and Treg), CD8+ T cells, monocytes, dendritic cells, and mast cells. Due to alternative splicing, CD366 exists as membrane-bound and soluble forms. Galectin-9 has been identified as the ligand for CD366. In humans, this receptor negatively regulates CD4+ T cells, influencing the secretion of some Th1- and Th17-related cytokines. CD366 has also been implicated in tolerance, autoimmune disease (e.g. multiple sclerosis), and HIV infection.

Applications Reported: This F38-2E2 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This F38-2E2 antibody has been pre-diluted and tested by flow cytometric analysis of stimulated normal human peripheral blood cells. This may be used at 5 µL (0.125 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test.


Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222-49) (100 µL of cell sample + 100 µL of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333-57) for up to 3 days in the dark at 4°C with minimal impact on

brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 488-561 nm; Emission: 667 nm; Laser: Blue Laser, Green Laser, Yellow-Green Laser

Product Images For CD366 (TIM3) Monoclonal Antibody (F38-2E2), PE-Cyanine5, eBioscience™

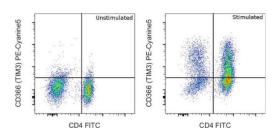

CD366 (TIM3) Antibody (15-3109-42) in Flow

Normal human peripheral blood cells were unstimulated (left) or stimulated for 72 hours with CD3 and CD28 Monoclonal Antibodies, Functional Grade (Product # 16-0037-85 and Product # 16-0289-85) (right). Cells were then stained with CD4 Monoclonal Antibody, FITC (Product # 11-0049-42) and CD366 (TIM3) Monoclonal Antibody, PE-Cyanine5. Total viable cells were used for analysis, as determined by Fixable Viability Dye eFluor 780 (Product # 65-0865-18).

View more figures on thermofisher.cn

□ 5 References

Western Blot (1)

The Biochemical journal

Phosphatidylserine binding directly regulates TIM-3 function.

"Published figure using CD366 (TIM3) monoclonal antibody (Product # 15-3109-42) in Western Blot" Authors: Smith CM,Li A,Krishnamurthy N,Lemmon MA

Year 2021

Flow Cytometry (4)

EBioMedicine

SLC1A1 mediated glutamine addiction and contributed to natural killer T-cell lymphoma progression with immunotherapeutic potential.

"Published figure using CD366 (TIM3) monoclonal antibody (Product # 15-3109-42) in Flow Cytometry"

Authors: Xiong J,Wang N,Zhong HJ,Cui BW,Cheng S,Sun R,Chen JY,Xu PP,Cai G,Wang L,Sun XJ,Huang JY,Zhao WL

Year 2021

Oncoimmunology

TIGIT* TIM-3* NK cells are correlated with NK cell exhaustion and disease progression in patients with hepatitis B virusrelated hepatocellular carcinoma.

"Published figure using CD366 (TIM3) monoclonal antibody (Product # 15-3109-42) in Flow Cytometry" Authors: Yu L,Liu X,Wang X,Yan F,Wang P,Jiang Y,Du J,Yang Z

Year 2021

View more Flow references on thermofisher.cn

More applications with references on thermofisher.cn

For Research Use Only, Not for use in diagnostic procedures. Not for resale without express authorization, Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CORPORMING PRODUCTS DURING THE WARRANTY PERFIGIO IS LIMITED. ACCIDENT OR FERIAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS, OLT STELLEY SOLE OPTION. THERE IS NO AGAINATION TO REPLAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS. AT SELLER'S SOLE OPTION. THERE IS NO AGAINATION TO REPLAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or vivy or in vivo therapeutic uses, or vivy or in vivo therapeutic uses, or vivy or in vivo therapeutic uses,