



CD73 Monoclonal Antibody (eBioTY/11.8 (TY/11.8)), Functional Grade, eBioscience™

Product Details	
Size	100 μg
Species Reactivity	Mouse
Published Species	Mouse
Host/Isotype	Rat / IgG1
Recommended Isotype Control	Rat IgG1 kappa Isotype Control (eBRG1), Functional Grade, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	eBioTY/11.8 (TY/11.8)
Conjugate	Functional Grade
Form	Liquid
Concentration	1 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	no preservative
Storage conditions	4° C
RRID	AB_1907408

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	4 Publications
Immunocytochemistry (ICC/IF)	-	2 Publications
Flow Cytometry (Flow)	0.5 µg/test	19 Publications
Functional Assay (FN)	Assay-Dependent	-
Control (Ctrl)	Assay-Dependent	-

Product Specific Information

Description: eBioTY/11.8 recognizes the with CD73 a 69-kDa GPI-anchored cell-surface protein with ecto-5'-nucleotidase activity. Expression on myeloid cells (CD11b) is restricted to the bone marrow. In human CD73 can be induced to secrete a soluble form with IL-2 suggesting a role in mediating activation signals. Differences between human and mouse CD73 have been reported. BALB/c mice have more CD4+CD73+ than CD8+CD73+ while the reciprocal is documented in humans.

CD73 is expressed on a subset of lymphocytes and increases during lymphocyte maturation. Recently it has been found that memory CD4 T cells express and are similar to the uncommitted primed precursor helper cells (Thpp) that can differentiate into Th1 or Th2 cells. Furthermore CD73 has been found on regulatory T cells.

Applications Reported: This eBioTY/11.8 (TY/11.8) antibody has been reported for use in flow cytometric analysis.

Applications Tested: This eBioTY/11.8 (TY/11.8) antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 0.5 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Storage and handling: Use in a sterile environment.


Filtration: 0.2 µm post-manufacturing filtered.

Purity: Greater than 90%, as determined by SDS-PAGE.

Endotoxin Level: Less than 0.001 ng/µg antibody, as determined by LAL assay.

Aggregation: Less than 10%, as determined by HPLC.

Product Images For CD73 Monoclonal Antibody (eBioTY/11.8 (TY/11.8)), Functional Grade, eBioscience™

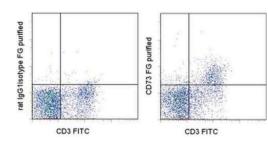

CD73 Antibody (16-0731-82) in Flow

Staining of BALB/c splenocytes with Anti-Mouse CD3e FITC (Product # 11-0031-82) and 0.25 µg of Rat IgG1 K Isotype Control Functional Grade Purified (Product # 16-4301-81) (left) or 0.25 µg of Anti-Mouse CD73 Functional Grade Purified (right) followed by Anti-Rat IgG PE (Product # 12-4822). Total viable cells were used for analysis.

View more figures on thermofisher.cn

□ 25 References

Immunohistochemistry (4)

Science advances

Dysregulation of ectonucleotidase-mediated extracellular adenosine during postmenopausal bone loss.

"Published figure using CD73 monoclonal antibody (Product # 16-0731-82) in Immunofluorescence"

Authors: Shih YV,Liu M,Kwon SK,lida M,Gong Y,Sangaj N,Varghese S

Year 2019

The Journal of experimental medicine

VEGF expands erythropoiesis via hypoxia-independent induction of erythropoietin in noncanonical perivascular stromal cells.

"Published figure using CD73 monoclonal antibody (Product # 16-0731-82) in Immunohistochemistry" Authors: Greenwald AC,Licht T,Kumar S,Oladipupo SS,Iyer S,Grunewald M,Keshet E

Year 2019

View more IHC references on thermofisher.cn

Immunocytochemistry (2)

The Journal of experimental medicine

VEGF expands erythropoiesis via hypoxia-independent induction of erythropoietin in noncanonical perivascular stromal cells.

"Published figure using CD73 monoclonal antibody (Product # 16-0731-82) in Immunohistochemistry" Authors: Greenwald AC,Licht T,Kumar S,Oladipupo SS,Iyer S,Grunewald M,Keshet E

Year 2019

Chemical senses

Characterization and turnover of CD73/IP(3)R3-positive microvillar cells in the adult mouse olfactory epithelium.

"16-0731 was used in Immunofluorescence to indicate that CD73-microvillar cells likely differentiate from proliferating progenitor cells and have a slow turnover."

Authors: Pfister S,Dietrich MG,Sidler C,Fritschy JM,Knuesel I,Elsaesser R

Year 2012

Species Mouse

More applications with references on thermofisher.cn

Flow (19)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation specifications and/or accompanying package inserts ("Documentation"), No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is used to properly trained individuals. Unless otherwise stated in the Documentation, this warranty pack are to extend to anyone other than the Buyer. Any mode of or sample increase in series illustrative of the general pend quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.
BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACE OR REFUND FOR FRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MALEJUER, (II) WISASTER OR EVENT OF FORCE MALEJUER, (II) WISASTER OR EVENT OF FORCE MALEJUER, (II) WISASTER OR EVENT OF FORCE MALEJUER (III) WISASTER OR EVENT OF FORCE MALEJUER, (III) WISASTER OR EVENT OF FORCE MALEJUER (III) WISASTER OR EVENT OF FORCE MALEJUENCE (III) OR SO THE PRODUCTS IN ANAMER FOR WHICH THEY WERE NOT DESIGNED, (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS IN ANAMER FOR WHICH THEY WERE NOT DESIGNED (IV) IMPROPER STORAGE AND HANDLING OF IN VIVO therapeutic uses, or any type of consumption by or application to human or animals.