Integrin alpha 4 beta 7 (LPAM-1) Monoclonal Antibody (DATK32 (DATK-32)), APC, eBioscience™

i i oudot Dotalio	
Size	100 µg
Species Reactivity	Mouse
Published Species	Mouse
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), APC, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	DATK32 (DATK-32)
Conjugate	APC
Excitation/Emission Max	651/660 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_1210577

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.25 µg/test	10 Publications

Product Specific Information

Description: The DATK32 monoclonal antibody binds to a combinatorial epitope of mouse LPAM-1 (alpha 4/ beta 7), which is a heterodimer of the the 154 kDa alpha 4 subunit and the 130 kDa beta 7 subunit. In mouse, the beta 7 subunit is found on the majority of mature lymphocytes and a small population of thymocytes and bone marrow cells. LPAM-1 is an integrin involved in the transendothelial migration of lymphocytes across high endothelial venules (HEV), and is know to bind several ligands including MAdCAM-1, VCAM-1 and fibronectin. Binding of the DATK32 monoclonal antibody has been shown to induce aggregation of the CD8+ T cell lymphoma TK1, and block LPAM-1-mediated cell adhesion.

Applications Reported: This DATK32 (DATK-32) antibody has been reported for use in flow cytometric analysis.


Applications Tested: This DATK32 (DATK-32) antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 0.25 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Excitation: 633-647 nm; Emission: 660 nm; Laser: Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.

1

Product Images For Integrin alpha 4 beta 7 (LPAM-1) Monoclonal Antibody (DATK32 (DATK-32)), APC, eBioscience™

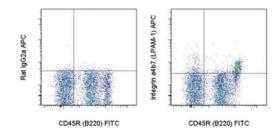

Integrin alpha 4 beta 7 (LPAM-1) Antibody (17-5887-82) in Flow Staining of C57Bl/6 bone marrow cells with Anti-Human/Mouse CD45R (B220) FITC (Product # 11-0452-82) and 0.125 µg of Rat IgG2a K Isotype Control APC (Product # 17-4321-81) (left) or 0.125 µg of Anti-Mouse Integrin alpha 4 beta 7 (LPAM-1) APC (right). Total viable cells were used for analysis.

□ 10 References

Flow Cytometry (10)

Immunity Epithelial colonization by gut dendritic cells promotes their functional diversification. "17-5887-82 was used in Flow cytometry/Cell sorting to investigate How tissue-specific cues orchestrate Dendritic cells diversification."	Year 2022 Species Mouse
Authors: Rivera CA,Randrian V,Richer W,Gerber-Ferder Y,Delgado MG,Chikina AS,Frede A,Sorini C,Maurin M, Kammoun-Chaari H,Parigi SM,Goudot C,Cabeza-Cabrerizo M,Baulande S,Lameiras S,Guermonprez P,Reis e Sousa C, Lecuit M,Moreau HD,Helft J,Vignjevic DM,Villablanca EJ,Lennon-Duménil AM	
Immunity An Id2 ^{RFP} -Reporter Mouse Redefines Innate Lymphoid Cell Precursor Potentials. "17-5887-82 was used in Flow Cytometry to suggest a revised model for ILC differentiation that redefines the cell-fate potential of helper-ILC-restricted Zbtb16+ ILCPs."	Year 2019 Species Mouse
Authors: Xu W,Cherrier DE,Chea S,Vosshenrich C,Serafini N,Petit M,Liu P,Golub R,Di Santo JP	

View more Flow references on thermofisher.cn

More applications with references on thermofisher.cn

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of said, as set forth in the Production documentation, specifications and effect at the time of said on substantially in conformance with published Product specifications in effect at the time of said set. as set forth in the Production documentation, this warranty specifications and effect at the time of said set. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not attend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR INVELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE ENterboy FOR NON-CONFINING PRODUCTS DURING THE WARRANTE NO CLEADING TO TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLED OFTION. THERE IS NO BULECATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLED OFTION CONFINING PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, [UI], USE OF THE PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, IC ON REGULERCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND in vitro diagnostic uses, or vitor in vito therapeutic uses, or any type of consumption by or application to human or animals.