



# Integrin alpha 4 beta 7 (LPAM-1) Monoclonal Antibody (DATK32 (DATK-32)), APC, eBioscience™

<b>Product Details</b>	
Size	100 μg
Species Reactivity	Mouse
Published Species	Mouse
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), APC, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	DATK32 (DATK-32)
Conjugate	APC
Excitation/Emission Max	651/660 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_1210577

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.25 µg/test	10 Publications

#### **Product Specific Information**

Description: The DATK32 monoclonal antibody binds to a combinatorial epitope of mouse LPAM-1 (alpha 4/ beta 7), which is a heterodimer of the the 154 kDa alpha 4 subunit and the 130 kDa beta 7 subunit. In mouse, the beta 7 subunit is found on the majority of mature lymphocytes and a small population of thymocytes and bone marrow cells. LPAM-1 is an integrin involved in the transendothelial migration of lymphocytes across high endothelial venules (HEV), and is know to bind several ligands including MAdCAM-1, VCAM-1 and fibronectin. Binding of the DATK32 monoclonal antibody has been shown to induce aggregation of the CD8+ T cell lymphoma TK1, and block LPAM-1-mediated cell adhesion.


Applications Reported: This DATK32 (DATK-32) antibody has been reported for use in flow cytometric analysis.


Applications Tested: This DATK32 (DATK-32) antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 0.25  $\mu$ g per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Excitation: 633-647 nm; Emission: 660 nm; Laser: Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.

## Product Images For Integrin alpha 4 beta 7 (LPAM-1) Monoclonal Antibody (DATK32 (DATK-32)), APC, eBioscience™

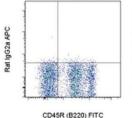

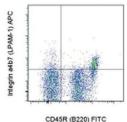

Integrin alpha 4 beta 7 (LPAM-1) Antibody (17-5887-82) in Flow Staining of C57Bl/6 bone marrow cells with Anti-Human/Mouse CD45R (B220) FITC (Product # 11-0452-82) and 0.125  $\mu g$  of Rat IgG2a K Isotype Control APC (Product # 17-4321-81) (left) or 0.125  $\mu g$  of Anti-Mouse Integrin alpha 4 beta 7 (LPAM-1) APC (right). Total viable cells were used for analysis.

### □ 10 References

### Flow Cytometry (10)

#### **Immunity**

## Epithelial colonization by gut dendritic cells promotes their functional diversification.

"17-5887-82 was used in Flow cytometry/Cell sorting to investigate How tissue-specific cues orchestrate Dendritic cells diversification."

Authors: Rivera CA,Randrian V,Richer W,Gerber-Ferder Y,Delgado MG,Chikina AS,Frede A,Sorini C,Maurin M, Kammoun-Chaari H,Parigi SM,Goudot C,Cabeza-Cabrerizo M,Baulande S,Lameiras S,Guermonprez P,Reis e Sousa C, Lecuit M,Moreau HD,Helft J,Vignjevic DM,Villablanca EJ,Lennon-Duménil AM

### Year

Species Mouse

#### Immunity

# An Id2<sup>RFP</sup>-Reporter Mouse Redefines Innate Lymphoid Cell Precursor Potentials.

"17-5887-82 was used in Flow Cytometry to suggest a revised model for ILC differentiation that redefines the cell-fate potential of helper-ILC-restricted Zbtb16+ ILCPs."

Authors: Xu W,Cherrier DE,Chea S,Vosshenrich C,Serafini N,Petit M,Liu P,Golub R,Di Santo JP

### **Year** 2019

Species Mouse

View more Flow references on thermofisher.cn

### More applications with references on thermofisher.cn

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization, Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production of specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, in this warranty provided herein is valid only when used by prospect pyrained individuals. Unless otherwise stated in the Documentation, in this warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED OR PREFUND FOR OR FERUND FOR THE NON-CONFORMING PRODUCTS, AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AT THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) WISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPES STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or vivy or in vivo therapeutic uses, or any type of consumption to hourand or animal tone to home and animal tone to the unauthorized commercial uses, in vitro diagnostic