# CD47 Monoclonal Antibody (B6H12), PerCP-eFluor™ 710, eBioscience™

Product Details	
Size	100 Tests
Species Reactivity	Human
Published Species	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), PerCP-eFluor™ 710, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	B6H12
Conjugate	PerCP-eFluor™ 710
Excitation/Emission Max	482/708 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_1944381

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.5 μg)/test	10 Publications

#### **Product Specific Information**

Description: The monoclonal antibody B6H12 reacts to CD47 also known as integrin-associated protein (IAP), and neurophilin. CD47 is a glycosylated five transmembrane protein with a small alternatively spliced cytoplasmic domain. CD47 is involved in adhesion through interactions with SIRP (signal regulator protein) and is non-covalently associated with beta3 integrins CD51 /CD61 and CD41/CD61. Furthermore this interaction can mediate bi-directional signaling to modify neural synaptic activity and regulate the phagocytic activities of macrophages. CD47 is the receptor for thrombospondin. T cell expression of CD47 can mediate activation or apoptosis (in the presence of high levels of thrombospondin). Recently stimulation of CD47 by monoclonal antibody has been shown to induce CD4+CD25- suppressive activity also increasing expression of Foxp3. Expression is found in the majority of hematopoietic cells including T and B cells, monocytes, platelets and erythrocytes (as part of the Rh complex). Expression is also found in non-hematopoietic cells.

This antibody has been reported to have neutralizing activity.

Applications Reported: This B6H12 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This B6H12 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5  $\mu$ L (0.5  $\mu$ g) per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells /test.

1


PerCP-eFluor® 710 emits at 710 nm and is excited with the blue laser (488 nm); it can be used in place of PerCP-Cyanine5.5. We recommend using a 710/50 bandpass filter, however, the 695/40 bandpass filter is an acceptable alternative. Please make sure that your instrument is capable of detecting this fluorochrome.

Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 µL cell sample + 100 µL IC Fixation Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency /compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 488 nm; Emission: 710 nm; Laser: Blue Laser.

Filtration: 0.2 µm post-manufacturing filtered.

## Product Images For CD47 Monoclonal Antibody (B6H12), PerCP-eFluor™ 710, eBioscience™

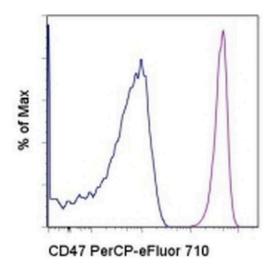

#### CD47 Antibody (46-0479-42) in Flow

Staining of normal human peripheral blood cells with Mouse IgG1 K Isotype Control PerCP-eFluor® 710 (Product # 46-4714-82) (blue histogram) or Anti-Human CD47 PerCP-eFluor® 710 (purple histogram). Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.cn

#### **□** 10 References

#### Flow Cytometry (10)

Pharmaceuticals (Basel, Switzerland)
Altered Membrane Expression and Function of CD11b Play a Role in the
Immunosuppressive Effects of Morphine on Macrophages at the
Nanomolar Level.

"Published figure using CD47 monoclonal antibody (Product # 46-0479-42) in Flow Cytometry" Authors: Yu PC,Hao CY,Fan YZ,Liu D,Qiao YF,Yao JB,Li CZ,Yu Y

# Oncoimmunology Year Schweinfurthin induces ICD without ER stress and caspase activation. 2022 "Published figure using CD47 monoclonal antibody (Product # 46-0479-42) in Flow Cytometry" Authors: Zhang R,Neighbors JD,Schell TD,Hohl RJ

View more Flow references on thermofisher.cn

**Year** 2023

### More applications with references on thermofisher.cn

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is provided herein is valid only when used by properly trained individuals. Unless otherwise of the general type and quality of goods and does not represent that any Product will conform to such model or sample timeled usage. This warranty to clean to textend to anyone other than the Buyer. Any model or sample timeles to subject to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample timeles to subject to normalize, EXPRESS OR INN INFRINGEMENT. BUYER'S EXCUSIVE TEMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTES, EXPRESS OR INN INFRINGEMENT. BUYER'S EXCUSIVE TEMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTE IN CEPIAIR. REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS) AT SELLER'S SOLE OPTION. THERE IS NO SBLEATION TO REPRESS OLE OPTION TO REPLAIR. REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS) AT SELLER'S SOLE OPTION. THERE IS NO SBLEATION TO REPRESS OLE OPTION PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLEGENCE OF OR BY BUYER', (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING