CD150 Monoclonal Antibody (mShad150), PerCP-eFluor™ 710, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Published Species	Mouse
Host/Isotype	Rat / IgG2b, lambda
Recommended Isotype Control	Rat IgG2b kappa Isotype Control (eB149/10H5), PerCP-eFluor™ 710, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	mShad150
Conjugate	PerCP-eFluor™ 710
Excitation/Emission Max	482/708 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2016699

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.125 µg/test	7 Publications

Product Specific Information

Description: The mShad150 monoclonal antibody reacts with mouse CD150, an ~70 kDa transmembrane glycoprotein also known as Signaling Lymphocyte Activation Molecule (SLAM). CD150 is expressed by T cells, in particular Th1, and B cells; this expression is rapidly upregulated upon activation. Immature thymocytes and dendritic cells also express this antigen. Signaling through SLAM in T cells induces proliferation and augmentation of the interferon-gamma response. Furthermore, SLAM is thought to play a role in adhesion between the T cell and antigen-presenting cell. The mShad150 antibody has been reported to also stain the CD150+CD244-CD48-CD41- population of pluripotent hematopoietic stem cells.

Applications Reported: This mShad150 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This mShad150 antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 0.125 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

PerCP-eFluor® 710 can be used in place of PE-Cy5, PE-Cy5.5 or PerCP-Cy5.5. PerCP-eFluor® 710 emits at 710 nm and is excited with the blue laser (488 nm). Please make sure that your instrument is capable of detecting this fluorochrome. For a filter configuration, we recommend using the 685 LP dichroic mirror and 710/40 band pass filter, however the 695/40 band pass filter is an acceptable alternative.


1

Our testing indicates that PerCP-eFluor® 710 conjugated antibodies are stable when stained samples are exposed to freshly prepared 2% formaldehyde overnight at 4°C, but please evaluate for alternative fixation protocols.

Excitation: 488 nm; Emission: 710 nm; Laser: Blue Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For CD150 Monoclonal Antibody (mShad150), PerCP-eFluor™ 710, eBioscience™

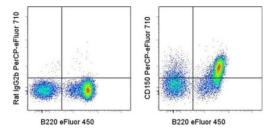

CD150 Antibody (46-1502-82) in Flow

Staining of C57BL/6 splenocytes with Anti-Human/Mouse CD45R (B220) eFluor® 450 (Product # 48-0452-82) and 0.06 µg of Rat IgG2b K Isotype Control PerCP-eFluor® 710 (Product # 46-4031-82) (left) or 0.06 µg of Anti-Mouse CD150 PerCP-eFluor® 710 (right). Total viable cells were used for analysis.

View more figures on thermofisher.cn

□ 7 References	
low Cytometry (7)	
Wellcome open research Activation of regulatory T cells triggers specific changes in glycosylation associated with Siglec-1-dependent inflammatory responses.	Year 2022
"Published figure using CD150 monoclonal antibody (Product # 46-1502-82) in Flow Cytometry"	
Authors: Wu G,Murugesan G,Nagala M,McCraw A,Haslam SM,Dell A,Crocker PR The Journal of experimental medicine	Year
Loss of tRNA-modifying enzyme Elp3 activates a p53-dependent	2021
antitumor checkpoint in hematopoiesis.	Species Mouse
"46-1502-82 was used in flow cytometry to investigate Elp3 and p53 in hematopoiesis."	MOUSE
Authors: Rosu A,EI Hachem N,Rapino F,Rouault-Pierre K,Jorssen J,Somja J,Ramery E,Thiry M,Nguyen L,Jacquemyn M,Daelemans D,Adams CM,Bonnet D,Chariot A,Close P,Bureau F,Desmet CJ	

View more Flow references on thermofisher.cn

More applications with references on thermofisher.cn

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided here in is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is minited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR INNELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, ITNESS FOR ANY PARTICULTS DURING THE WARRANT PERIOD IS LIMITED TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAR, REPLACE COR REFUND FOR THE NON-CONFORMING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION. THERE IS NO BULGATION TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION. THERE IS NO BULGATION TO REPARR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation, accompanying the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorizad commercial uses, in vitro diagnostic uses, or vivo or in vivo therapueutic uses, or any type of consumption to human or animals.

2