B7-H6 Monoclonal Antibody (JAM1EW), PerCP-eFluor[™] 710, eBioscience[™]

Product Details	
Size	100 Tests
Species Reactivity	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), PerCP-eFluor™ 710, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	JAM1EW
Conjugate	PerCP-eFluor™ 710
Excitation/Emission Max	482/708 nm
Immunogen	B7-H6 recombinant protein
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2784692

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.5 μg)/test	-

Product Specific Information

Description: This JAM1EW monoclonal antibody recognizes human B7-H6 (NCR3LG1) protein on the surface of cells. It recognizes protein on live but not on paraformaldehyde-fixed cells. The JAM1EW clone is specific to B7-H6. When tested in flow cytometry against cells transfected with either B7-H4, B7-H5, B7-H6, or B7-H7, clone JAM1EW only stains positively on B7-H6 transfected cells.

Applications Tested: This JAM1EW antibody has been pre-diluted and tested by flow cytometric analysis of the HL-60 cell line. This may be used at 5 μ L (0.5 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.

PerCP-eFluor 710 emits at 710 nm and is excited with the blue laser (488 nm); it can be used in place of PerCP-Cyanine5.5. We recommend using a 710/50 bandpass filter, however, the 695/40 bandpass filter is an acceptable alternative. Please make sure that your instrument is capable of detecting this fluorochrome.

Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.


Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222-49) (100 μ L of cell sample + 100 μ L of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333-57) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be

1

made, but clone specific performance should be determined empirically.


Excitation: 488 nm; Emission: 710 nm; Laser: Blue Laser

Product Images For B7-H6 Monoclonal Antibody (JAM1EW), PerCP-eFluor™ 710, eBioscience™

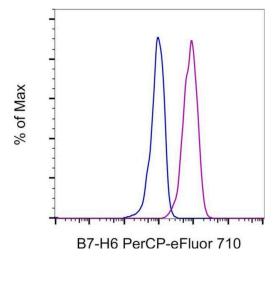

B7-H6 Antibody (46-6526-42) in Flow

HeLa cells were stained with Mouse IgG1 kappa Isotype Control, PerCP-eFluor 710 (Product # 46-4714-82) (blue) or B7-H6 Monoclonal Antibody, PerCP-eFluor 710 (purple). Total viable cells, as determined by Fixable Viability Dye eFluor 450 (Product # 65-0863-18), were used for analysis.

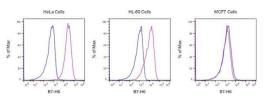

B7-H6 Antibody (46-6526-42)

Surface staining of cell lines. As expected based on known relative expression patterns, B7-H6 clone JAM1EW stains HeLa and HL-60, but not MCF7 cell lines. Details: HeLa (left histogram), HL-60 (middle histogram), or MCF7 (right histogram) cells were stained with Mouse IgG1 kappa Isotype Control (blue) (clone P3.6.2.8.1) or B7-H6 Monoclonal Antibody (clone JAM1EW, purple). Total viable cells, as determined by Fixable Viability Dye eFluor 450 (Product # 65-0863-18), were used for analysis. {RE}