Performance guarenteed

CD38 Monoclonal Antibody (HIT2), APC-eFluor™ 780, eBioscience™

Product Details	
Size	100 Tests
Species Reactivity	Human
Published Species	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), APC-eFluor™ 780, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	HIT2
Conjugate	APC-eFluor™ 780
Excitation/Emission Max	756/785 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_11218674

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.25 μg)/test	12 Publications

Product Specific Information

Description: The HIT2 monoclonal antibody reacts with the human CD38 molecule, an approximately 45 kDa type II transmembrane protein. CD38 is expressed by thymocytes, peripheral B cells including plasma cells, activated T cells, and monocytes. CD38 is a counter-receptor for CD31 and has both cyclase and hydrolase enzymatic activity.

Applications Reported: This HIT2 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This HIT2 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 μ L (0.25 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.

APC-eFluor 780 emits at 780 nm and is excited with the Red laser (633 nm). Please make sure that your instrument is capable of detecting this fluorochome.

Light sensitivity: This tandem is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.


Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 μ L cell sample + 100 μ L IC Fixation Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency /compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific

performance should be determined empirically.

Excitation: 633-647 nm; Emission: 780 nm; Laser: Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For CD38 Monoclonal Antibody (HIT2), APC-eFluor™ 780, eBioscience™

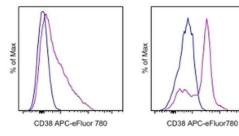

CD38 Antibody (47-0389-42) in Flow

Staining of normal human peripheral blood cells with Mouse IgG1 K Isotype Control APC-eFluor® 780 (Product # 47-4714-82) (blue histogram) or Anti-Human CD38 APC-eFluor® 780 (purple histogram). Cells in the lymphocyte (left) or monocyte (right) gate were used for analysis.

View more figures on thermofisher.cn

A 12 References			
Year 2023			
Year			
2022			

View more Flow references on thermofisher.cn

More applications with references on thermofisher.cn

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. No OTHER WARRANTES, EXERGES OR IMPLED, ARE CRANTED INCLUDING WITHOUTUILIMITATION. IMPLEID WARRANTES OF MERCHANTABILITY, FITNESS FOR ANTE INCES, CP NON NHORENGE, CP RON NHORENGE, CP NON NHORENGE, CP RON N