Performanc guarenteed

TER-119 Monoclonal Antibody (TER-119), eFluor™ 450, eBioscience™

Product Details		
Size	100 µg	
Species Reactivity	Mouse	
Published Species	Fish, Mouse, Human	
Host/Isotype	Rat / IgG2b, kappa	
Recommended Isotype Control	Rat IgG2b kappa Isotype Control (eB149/10H5), eFluor™ 450, eBioscience™	
Class	Monoclonal	
Туре	Antibody	
Clone	TER-119	
Conjugate	eFluor™ 450	
Excitation/Emission Max	405/445 nm	
Form	Liquid	
Concentration	0.2 mg/mL	
Purification	Affinity chromatography	
Storage buffer	PBS, pH 7.2	
Contains	0.09% sodium azide	
Storage conditions	4° C, store in dark, DO NOT FREEZE!	
RRID	AB_1518808	

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	6 Publications
Immunocytochemistry (ICC/IF)	-	1 Publication
Flow Cytometry (Flow)	0.5 µg/test	77 Publications

Product Specific Information

Description: The TER-119 monoclonal antibody reacts with mouse erythroid cells from early proerythroblast to mature erythrocyte stages. The TER-119 antigen is present in yolk sac, fetal and newborn liver, but is not expressed by cells carrying BFU-E and CFU-E activities. Several erythroleukemia cell lines tested so far are negative for expression of TER-119 antigen even after dimethylsulfoxide stimulation. Biochemical and molecular analysis of the TER-119 antigen indicate that this molecule is associated with the surface glycophorin A, but is not a typical glycophorin.

Applications Reported: This TER-119 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This TER-119 antibody has been tested by flow cytometric analysis of mouse bone marrow cells. This can be used at less than or equal to 0.5 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.


eFluor® 450 is an alternative to Pacific Blue®. eFluor® 450 emits at 445 nm and is excited with the Violet laser (405 nm). Please make sure that your instrument is capable of detecting this fluorochome.

1

Excitation: 405 nm; Emission: 445 nm; Laser: Violet Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For TER-119 Monoclonal Antibody (TER-119), eFluor™ 450, eBioscience™

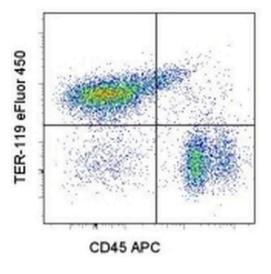

TER-119 Antibody (48-5921-82) in Flow

Staining of C57BL/6 bone marrow cells with Anti-Mouse CD45 APC (Product # 17-0451-82) and 0.25 μ g of Anti-Mouse TER-119 eFluor® 450. Total viable cells were used for analysis.

View more figures on thermofisher.cn

□84 References

Immunohistochemistry (6)

The Journal of clinical investigation Ribonuclease inhibitor 1 regulates erythropoiesis by controlling GATA1 translation.

"Published figure using TER-119 monoclonal antibody (Product # 48-5921-82) in Immunocytochemistry"

Authors: Chennupati V, Veiga DF, Maslowski KM, Andina N, Tardivel A, Yu EC, Stilinovic M, Simillion C, Duchosal MA, Quadroni M,Roberts I,Sankaran VG,MacDonald HR,Fasel N,Angelillo-Scherrer A,Schneider P,Hoang T,Allam R

Development (Cambridge, England)

Transient loss of venous integrity during developmental vascular remodeling leads to red blood cell extravasation and clearance by lymphatic vessels.

"Published figure using TER-119 monoclonal antibody (Product # 48-5921-82) in Immunofluorescence" Authors: Zhang Y, Daubel N, Stritt S, Mäkinen T

Immunocytochemistry (1)

The Journal of clinical investigation Year 2018 Ribonuclease inhibitor 1 regulates erythropoiesis by controlling GATA1 translation. "Published figure using TER-119 monoclonal antibody (Product # 48-5921-82) in Immunocytochemistry" Authors: Chennupati V, Veiga DF, Maslowski KM, Andina N, Tardivel A, Yu EC, Stilinovic M, Simillion C, Duchosal MA,

Quadroni M,Roberts I,Sankaran VG,MacDonald HR,Fasel N,Angelillo-Scherrer A,Schneider P,Hoang T,Allam R

Flow Cytometry (77)

PloS one

Engineered red blood cells carrying PCSK9 inhibitors persistently lower LDL and prevent obesity.

"Published figure using TER-119 monoclonal antibody (Product # 48-5921-82) in Flow Cytometry" Authors: Deshycka R,Sudaryo V,Huang NJ,Xie Y,Smeding LY,Choi MK,Ploegh HL,Lodish HF,Pishesha N

View more Flow references on thermofisher cn

More applications with references on thermofisher.cn

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation For Research Use Only. Not for use in alignostic procedures. Not for resiste without express authorization. Products are warranted to operate or perform substantiation in conformance with published Product specifications in effect at the time of salability for a substantiation. If a contract are marranted to appear to perform substantiation in the product specifications and the set of salability for a substantiation. If a contract are performed and the accumentation, and of accompanying package inserts ("Documentation," ho claim of suitability for use in applications is meter ("Based The warranty for added herein is valid only when used by properly trainations in meter ("Based The warranty for added herein is valid only when used by properly trainations in the Documentation, his warranty for added the reline is valid only when used by properly trainations in meter ("Based The warranty for added herein is valid only when used by properly trainations in the Documentation, his warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. No OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, ON NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEMPTOR NON CONTES OUTS DURING THE WARRANTY EXPRISION IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE RON-CONFORMING PRODUCTS AS THE RESULT OF (") ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (") IMPROPER STORAGE AND MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (") IMPROPER STORAGE AND MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (") IMPROPER STORAGE AND MANNER FOR WHICH THEY WERE NOT DESIGNED, or (v) IMPROPER STORAGE AND MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (v) IMPROPER STORAGE AND MANNER FOR WHIC HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the docum tation accompanying the Product, the Product is int ed for research only and is not to be us ed for any other purpose, including without limitati diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals

Year 2018

Year 2018

View more IHC references on thermofisher.cn

Year 2021