# CD41a Monoclonal Antibody (HIP8), PE-eFluor™ 610, eBioscience™

Product Details	
Size	100 Tests
Species Reactivity	Human
Published Species	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), PE-eFluor™ 610, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	HIP8
Conjugate	PE-eFluor™ 610
Excitation/Emission Max	565/606 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2848488

Applications	Tested Dilution	Publications
Western Blot (WB)	-	1 Publication
Flow Cytometry (Flow)	5 μL (0.06 μg)/test	3 Publications

#### **Product Specific Information**

Description: The HIP8 monoclonal antibody reacts with the human CD41 molecule, the integrin alpha IIb also known as platelet GPIIb. CD41 non-covalently associates with integrin beta3 (GPIIIa, CD61) and is expressed by megakaryocytes and platelets. The CD41/CD61 complex is a receptor for fibronectin, fibrinogen, von Willebrand factor, vitronectin and thrombospondin and mediates platelets aggregation. HIP8 blocks platelet aggregation.

Applications Reported: This HIP8 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This HIP8 antibody has been pre-diluted and tested by flow cytometric analysis of human platelets. This may be used at 5  $\mu$ L (0.06  $\mu$ g) per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.

PE-eFluor 610 can be excited with laser lines from 488-561 nm and emits at 607 nm. We recommend using a 610/20 band pass filter (equivalent to PE-Texas Red). Please make sure that your instrument is capable of detecting this fluorochrome.


Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.

1

Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222-49) (100 µL of cell sample + 100 µL of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333-57) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 488-561 nm; Emission: 607 nm; Laser: Blue Laser, Green Laser, Yellow-Green Laser

### Product Images For CD41a Monoclonal Antibody (HIP8), PE-eFluor™ 610, eBioscience™

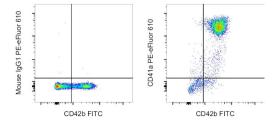

#### CD41a Antibody (61-0419-42) in Flow

Human platelets were stained with CD42b Monoclonal Antibody, FITC (Product # 11-0429-42) and Mouse IgG1 kappa Isotype Control, PE-eFluor 610 (Product # 61-4714-82) (left) or CD41a Monoclonal Antibody, PE-eFluor 610 (right). Total cells were used for analysis.

#### View more figures on thermofisher.cn

#### 4 References

#### Western Blot (1)

Communications biology	Year
A substrate-trapping strategy to find E3 ubiquitin ligase substrates	2020
identifies Parkin and TRIM28 targets.	Species
"61-0419-42 was used in Western Blot to develop an E3 ligase substrate-trapping technique that can identify substrates of E3 ligases."	Human
Authors: Watanabe M,Saeki Y,Takahashi H,Ohtake F,Yoshida Y,Kasuga Y,Kondo T,Yaguchi H,Suzuki M,Ishida H, Tanaka K,Hatakeyama S	1:2000

#### Flow Cytometry (3)

The Journal of experimental medicine Heterogeneous disease-propagating stem cells in juvenile myelomonocytic leukemia.		<b>Year</b> 2021	
"Published figure using CD41a monoclonal antibody (Product # 61-0419-42) in Flow Cytometry"			
	Authors: Louka E.Povinelli B.Rodriguez-Meira A.Buck G.Wen WX.Wang G.Sousos N.Ashlev N.Hamblin A.Booth CAG.		

Roy A,Elliott N,Iskander D,de la Fuente J,Fordham N,O'Byrne S,Inglott S,Norfo R,Salio M,Thongjuea S,Rao A,Roberts I, Mead AJ

#### Cancers

## Platelet-Derived GARP Induces Peripheral Regulatory T Cells-Potential Impact on T Cell Suppression in Patients with Melanoma-Associated Thrombocytosis.

"Published figure using CD41a monoclonal antibody (Product # 61-0419-42) in Flow Cytometry"

Authors: Zimmer N,Krebs FK,Zimmer S,Mitzel-Rink H,Kumm EJ,Jurk K,Grabbe S,Loquai C,Tuettenberg A

View more Flow references on thermofisher.cn

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No caim of subjected to normal, proper and intended usage. This warranty get on other than the Buyer. Any model or sample furnished to buyers buyer obtained to the property trained individuals. Unless otherwise stated in the Documentation, this warranty get from date of simple twinnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample furnished to Buyer is morely illustrative of the general type and quality of goods and does not erpresent that any Product will conform to such model or sample furnished to Buyer is morely illustrative of the general type and quality of goods and does not erpresent that any Product site ("DoroRNING PRODUCTS DURING THE WARRANTES, EXPRESS OR INVELD, VER CENT OF FOR ON-CONFORMING PRODUCTS DURING THE WARRANTES, EXPRESS OR INVELD, VER PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLEND THE WARRANTES NO EVENT OF FORCE MAJEURCE OR REFUND FOR REFUND FOR THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURCE, ION REGULARCE OR OR BY EXPRESS OR BY EXPRESS OR INVELON ON REGULARCE OR OR BY EXPRESS OR INVELON ON REGULARCE OR REFUND FOR THE RESOLUT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURCE, ION REGULARCE OR OR BY EVENT OF FORCE MAJEURCE, ION REGULARCE OR OR BY EVENT OF FORCE MAJEURCE, ION REGULARCE OR BY EVENT O

2

Year 2020