CD223 (LAG-3) Monoclonal Antibody (3DS223H), Super Bright[™] 436, eBioscience[™]

Product Details	
Size	100 Tests
Species Reactivity	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), Super Bright [™] 436, eBioscience [™]
Class	Monoclonal
Туре	Antibody
Clone	3DS223H
Conjugate	Super Bright™ 436
Excitation/Emission Max	413/431 nm
Immunogen	hCD223-FC PFA 4%
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2637155

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.125 μg)/test	7 Publications

Product Specific Information

Description: This 3DS223H monoclonal antibody recognizes human CD223 also known as Lymphocyte Activation Gene 3 (LAG-3). LAG-3 is a 70-kDa surface glycoprotein belonging to the Ig superfamily with homology to CD4. LAG-3 binds to MHC class II with higher affinity than CD4 and is thought to be involved in the negative regulation of T cell activation and homeostatic proliferation. Surface expression of LAG-3 has been reported on activated T cells (including regulatory T cells) and NK cells. CD8+ T cells usually express LAG-3 at significantly higher levels than CD4+ T cells. Coexpression of LAG-3 and CD49b has been proposed to identify human and mouse Type 1 regulatory T cells (Tr1 cells).

This 3DS223H antibody will recognize a formaldehyde-fixed epitope.

Applications Reported: This 3DS223H antibody has been reported for use in flow cytometric analysis.

Applications Tested: This 3DS223H antibody has been pre-titrated and tested by flow cytometric analysis of stimulated normal human peripheral blood cells. This can be used at 5 μ L (0.125 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.


Super Bright 436 can be excited with the violet laser line (405 nm) and emits at 436 nm. We recommend using a 450/50 bandpass filter, or equivalent. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.

Excitation: 405 nm; Emission: 436 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD223 (LAG-3) Monoclonal Antibody (3DS223H), Super Bright™ 436, eBioscience™

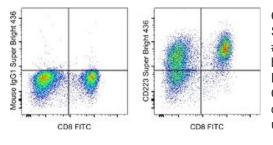

CD223 (LAG-3) Antibody (62-2239-42) in Flow

Staining of 3-day Anti-Human CD3 and CD28 Functional Grade Purified (Product # 16-0037-81 and Product # 16-0289-81)-stimulated normal human peripheral blood cells with Anti-Human CD8a FITC (Product # 11-0087-42) and Mouse IgG1 K Isotype Control Super Bright 436 (Product # 62-4714-82) (left) or Anti-Human CD223 (LAG-3) Super Bright 436 (right). Viable cells in the lymphocyte gate, as determined by 7-AAD Viability Staining Solution (Product # 00-6993-50) were used for analysis.

View more figures on thermofisher.cn

□ 7 References Flow Cytometry (7)			
Frontiers in immunology Bispecific T-Cell Engaging Antibodies Against MUC16 Demonstrate Efficacy Against Ovarian Cancer in Monotherapy and in Combination With PD-1 and VEGF Inhibition.	Year 2021		
"Published figure using CD223 (LAG-3) monoclonal antibody (Product # 62-2239-42) in Flow Cytometry" Authors: Yeku OO,Rao TD,Laster I,Kononenko A,Purdon TJ,Wang P,Cui Z,Liu H,Brentjens RJ,Spriggs D			

View more Flow references on thermofisher.cn

More applications with references on thermofisher.cn

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production occumentation, specifications and/or sample furning package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Indees otherwise stated in the Documentation, this warranty is united to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR INMELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, FITNESS FOR ANY PARTICURS DEVENDES, OR NON INFRINCEMENT. BUYER'S SCLUSVER EMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS OF NON-CONFORMING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS OF SOLE OFTION TO FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorizad commercial uses, in vitro diagnostic uses, or any type of consumption to human or animals.

2