

CD278 (ICOS) Monoclonal Antibody (ISA-3), Super Bright™ 436, eBioscience™

Product Details	
Size	100 Tests
Species Reactivity	Human
Published Species	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), Super Bright™ 436, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	ISA-3
Conjugate	Super Bright™ 436
Excitation/Emission Max	413/431 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2637393

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 µL (0.25 µg)/test	10 Publications
Neutralization (Neu)	-	1 Publication
T-Cell Activation (TCA)	-	2 Publications

Product Specific Information

Description: The ISA-3 monoclonal antibody reacts with human CD278 also known as Inducible Costimulatory molecule (ICOS), H4, CRP-1 and AILIM. ICOS is a T cell specific activation molecule and a third member of the CD28/CTLA-4 family. Human ICOS has a relative molecular mass of 55-60 kDa, composed of 27 kDa and 29 kDa chains. Human ICOS on activated T cells has potent costimulatory activity for T cell activation and is required for humoral immune responses, in particular for memory B cell and plasma cell generation. ICOS binds to its ligand, B7h/B7RP-1 expressed on activated APCs (antigen presenting cells) and on a number of inflamed peripheral tissues. Plate-bound ISA-3 is costimulatory for T cells and induces production of IL-4, IL-5, IL-10 and other cytokines, but not IL-2. ISA-3 has the same reactivity pattern and characteristics as F44. ISA-3 was generated against the human ICOS antigen. C398.4A, anti-mouse ICOS/H4 (Product # 14-9949-82), was shown to cross-react with human ICOS but binds to an epitope different from ISA-3. C398.4A stains activated cells brighter than ISA-3; however, it also exhibits higher staining of non-activated human peripheral blood or isolated PBMC. To achieve the brightest staining of ICOS on activated human T cells, please use Product # 13-9948-82 or Product # 12-9948-42 rather than Product # 11-9948-42 .

Applications Reported: This ISA-3 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This ISA-3 antibody has been pre-titrated and tested by flow cytometric analysis of stimulated normal

human peripheral blood cells. This can be used at 5 μ L (0.25 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.


Super Bright 436 can be excited with the violet laser line (405 nm) and emits at 436 nm. We recommend using a 450/50 bandpass filter, or equivalent. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.

Excitation: 405 nm; Emission: 436 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD278 (ICOS) Monoclonal Antibody (ISA-3), Super Bright™ 436, eBioscience™

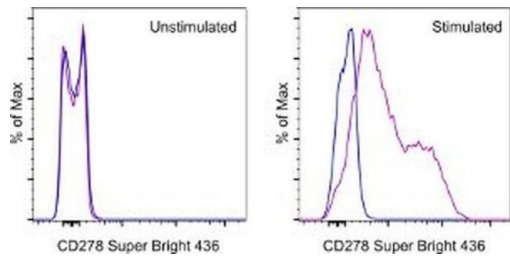

CD278 (ICOS) Antibody (62-9948-42) in Flow
Staining of unstimulated (left) or Anti-Human CD3 and CD28 Functional Grade Purified (Product # 16-0037-81 and Product # 16-0289-81)-stimulated (right) normal human peripheral blood cells with Mouse IgG1 K Isotype Control Super Bright 436 (Product # 62-4714-82) (blue histogram) or Anti-Human CD278 (ICOS) Super Bright 436 (purple histogram). Total viable cells, as determined by 7-AAD Viability Staining Solution (Product # 00-6993-50), were used for analysis.

[View more figures on thermofisher.cn](https://thermofisher.cn)

Flow Cytometry (10)

eLife	Year 2021
Impaired HA-specific T follicular helper cell and antibody responses to influenza vaccination are linked to inflammation in humans.	
"Published figure using CD278 (ICOS) monoclonal antibody (Product # 62-9948-42) in Flow Cytometry"	
Authors: Hill DL,Whyte CE,Innocentin S,Lee JL,Dooley J,Wang J,James EA,Lee JC,Kwok WW,Zand MS,Liston A,Carr EJ,Linterman MA	
Scientific reports	
Changes in expression of PD-L1 on peripheral T cells in patients with melanoma and lung cancer treated with PD-1 inhibitors.	
"Published figure using CD278 (ICOS) monoclonal antibody (Product # 62-9948-42) in Flow Cytometry"	
Authors: Dart SJ,Cook AM,Millward MJ,McDonnell AM,Chin WL,Hakeem MU,Meniawy TM,Bowyer SE	

[View more Flow references on thermofisher.cn](#)

Neutralization (1)

International journal of chronic obstructive pulmonary disease	Year 2022
Inducible Costimulator-C-X-C Motif Chemokine Receptor 3 Signaling is Involved in Chronic Obstructive Pulmonary Disease Pathogenesis.	
"Published figure using CD278 (ICOS) monoclonal antibody (Product # 62-9948-42) in Neutralization"	
Authors: Li DY,Chen L,Miao SY,Zhou M,Wu JH,Sun SW,Liu LL,Qi C,Xiong XZ	

T-Cell Activation (2)

Oncoimmunology	Year 2016
Polyfunctional Melan-A-specific tumor-reactive CD8(+) T cells elicited by dacarbazine treatment before peptide-vaccination depends on AKT activation sustained by ICOS.	
"Published figure using CD278 (ICOS) monoclonal antibody (Product # 62-9948-42) in T-Cell Activation"	
Authors: Franzese O,Palermo B,Di Donna C,Sperduti I,Ferraresi V,Stabile H,Gismondi A,Santoni A,Nisticò P	

More applications with references on thermofisher.cn

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.