CD62L (L-Selectin) Monoclonal Antibody (MEL-14), NovaFluor™ Yellow 660, eBioscience™

Product Details

Size	25 µg
Species Reactivity	Mouse
Host/Isotype	Rat / IgG2a, kappa
Class	Monoclonal
Туре	Antibody
Clone	MEL-14
Conjugate	NovaFluor™ Yellow 660
Excitation/Emission Max	550/664 nm
Form	Liquid
Concentration	0.1 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.2 µg/test	-

Product Specific Information

Description: The MEL-14 monoclonal antibody reacts with mouse CD62L, a 76 kDa member of the selectin family. CD62L is expressed by neutrophils, monocytes, and subsets of T, B, and NK cells and binds a number of glycosylated, fucosylated, sulfated sialylated glycoproteins including CD34, glycam-1 and MAdCam-1. These interactions mediate rolling of lymphocytes on activated endothelium at the sites of inflammation and homing of cells to the high endothelial venules (HEV) of peripheral lymphoid tissues.

Each product contains 1 vial of NovaFluor conjugate and 1 vial of CellBlox Plus Blocking Buffer .

Applications Reported: The MEL-14 antibody has been reported for use in flow cytometric analysis.

Applications Tested: The MEL-14 antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 0.2 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

NovaFluor dyes are not compatible with DNA intercalating viability dyes. Do not use viability dyes such as propidium iodide, 7actinomycin D (7-AAD) and DAPI. Invitrogen LIVE/DEAD Fixable Dead Cell stains are recommended for use with NovaFluor dyes.


This NovaFluor conjugate has been updated to ship with CellBlox Plus Blocking Buffer (Cat. No. (C001T06F01)). This buffer contains formulation improvements over CellBlox. CellBlox Plus Blocking Buffer is required for optimal staining with NovaFluor conjugates and should be used in all experiments where NovaFluor conjugates are used. Whenever possible, we recommend adding CellBlox Plus Blocking Buffer to antibody cocktails/master mixes prior to combining with cells. Add 5 μ L per sample (regardless of the number of NovaFluors in your panel) to use the antibody cocktail as intended. For single-color controls, use 5 μ L of CellBlox Blocking Buffer per 100 μ L of cell sample containing 10^3 to 10^8 cells.

NovaFluor conjugates are based on Phiton[™] technology utilizing novel nucleic acid dye structures that allow for engineered fluorescent signatures with consideration for spillover and spread impacts. Learn more

Excitation: 552 nm; Emission: 663 nm; Laser: 561 nm (Yellow) Laser


Our internal testing shows that NovaFluor Yellow 660 non-specifically stains B cells in SJL mice. Non-specific staining has not been observed in BALB/c or C57BL/6 mice. Other strains have not been tested. See the Antibody Testing Data for an example of this strain-dependent difference.

Product Images For CD62L (L-Selectin) Monoclonal Antibody (MEL-14), NovaFluor™ Yellow 660, eBioscience™

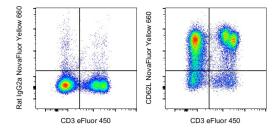

CD62L (L-Selectin) Antibody (M006T02Y04-A) in Flow


C57BL/6 mouse splenocytes were stained with CD3e Monoclonal Antibody, eFluor 450 (Product # 48-0031-82) and 0.2 µg of Rat IgG2a, kappa Isotype Control, NovaFluor Yellow 660 (left) or 0.2 µg of CD62L Monoclonal Antibody, NovaFluor Yellow 660 (right). Total viable cells were used for analysis, as determined by LIVE/DEAD Blue (Product # L34962).

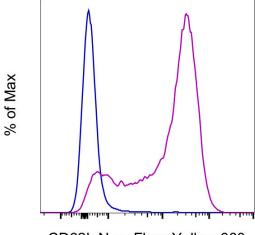

CD62L (L-Selectin) Antibody (M006T02Y04-A) in Flow

C57BL/6 mouse splenocytes were stained with 0.2 µg of Rat IgG2a, kappa Isotype Control, NovaFluor Yellow 660 (blue histogram) or 0.2 µg of CD62L Monoclonal Antibody, NovaFluor Yellow 660 (purple histogram) and acquired in the YG4 channel on a 5-laser Cytek Aurora. Cells in the splenocyte gate were used in the analysis.

CD62L NovaFluor Yellow 660

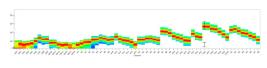

CD62L (L-Selectin) Antibody (M006T02Y04-A) in Flow

Spectral signature for NovaFluor Yellow 660 collected on a 5-laser Cytek Aurora Full Spectrum flow cytometer using Cytek assay settings. Human peripheral blood mononuclear cells were stained with anti-human CD4 (SK3) and signatures displayed following gating on the lymphocyte population.

View more figures on thermofisher.cn

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individus. Unsess otherwise stated in the Documentation, this warranty is imited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample turnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR INNELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEUP FOR NON-CONFORMING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS ON TO REPERING REPLAY. REPLACE OR REFUND FOR THE RON-CONFORMING PRODUCTS SOLE OPTION. THERE IS NO OBLICATION TO REPAR, REPLACE OR REFUND FOR THE RODUCTS AS THE RESULT OF (I) ACCLIEDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUBE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS UNES otherwise expressly stated on the Product in the documentation or animals.

2