

CD45 Polyclonal Antibody

Product Details	
Size	100 μL
Species Reactivity	Human, Mouse, Rat
Published Species	Human
Host/Isotype	Rabbit / IgG
Class	Polyclonal
Туре	Antibody
Conjugate	Unconjugated
Immunogen	Recombinant fusion protein containing a sequence corresponding to amino acids 194-390 of human CD45 (NP_002829.3)
Form	Liquid
Concentration	0.44 mg/mL
Purification	Affinity Chromatography
Storage buffer	PBS, pH 7.3, with 50% glycerol
Contains	0.05% ProClin 300
Storage conditions	-20° C, Avoid Freeze/Thaw Cycles
RRID	AB_2807863


Applications	Tested Dilution	Publications
Western Blot (WB)	1:500-1:1,000	-
Immunohistochemistry (Paraffin) (IHC (P))	1:50-1:200	-
Immunocytochemistry (ICC/IF)	1:50-1:200	-
Flow Cytometry (Flow)	-	1 Publication

Product Specific Information

Positive Samples: Mouse thymus; Cellular Location: Membrane, Membrane raft, Single-pass type I membrane protein


Immunogen sequence: SDAYLNASET TTLSPSGSAV ISTTTIATTP SKPTCDEKYA NITVDYLYNK ETKLFTAKLN VNENVECGNN TCTNNEVHNL TECKNASVSI SHNSCTAPDK TLILDVPPGV EKFQLHDCTQ VEKADTTICL KWKNIETFTC DTQNITYRFQ CGNMIFDNKE IKLENLEPEH EYKCDSEILY NNHKFTNASK IIKTDFG

Product Images For CD45 Polyclonal Antibody

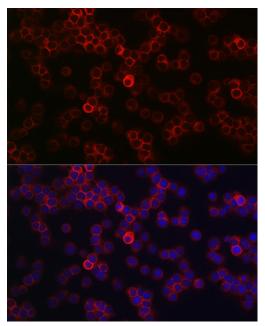


CD45 Antibody (PA5-96061) in ICC/IF

Immunofluorescence analysis of CD45 in Jurkat cells. Samples were incubated with CD45 Polyclonal antibody (Product # PA5-96061) using a dilution of 1:100 (40x lens). Blue: DAPI for nuclear staining.

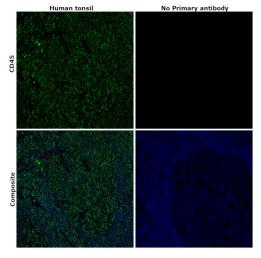

CD45 Antibody (PA5-96061) in IHC (P)

Immunohistochemical analysis of CD45 was performed using formalin-fixed paraffin-embedded human tonsil tissue sections. To expose the target protein, heat-induced epitope retrieval was performed on de-paraffinized sections using eBioscience™ IHC Antigen Retrieval Solution - High pH (10X) (Product # 00-4956-58) diluted to 1X solution in water in a microwave at 100 degree Celsius for 10 minutes. Following antigen retrieval, the sections were blocked with 2% normal goat serum in 1X PBS for 45 minutes at room temperature and then probed with or without CD45 Polyclonal Antibody (Product # PA5-96061) at 1: 100 dilution in 0.1% normal goat serum overnight at 4 degree Celsius in a humidified chamber. Detection was performed using Goat anti-Rabbit IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor Plus 488 (Product # A32731) at a dilution of 1:2000 in 0.1% normal goat serum for 45 minutes at room temperature. ReadyProbes™ Tissue Autofluorescence Quenching Kit (Product # R37630) was used to quench autofluorescence from the tissues. Nuclei were stained with DAPI (Product # D1306) and the sections were mounted using ProLong™ Glass Antifade Mountant (Product # P36984). The images were captured on EVOS™ M7000 Imaging System (Product # AMF7000) at 20X magnification and deconvoluted externally.

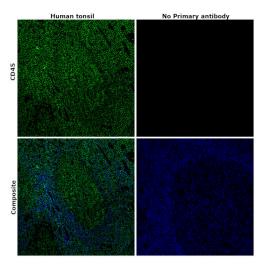

CD45 Antibody (PA5-96061) in IHC (P)

Immunohistochemical analysis of CD45 was performed using formalin-fixed paraffin-embedded human tonsil tissue sections. To expose the target protein, heat-induced epitope retrieval was performed on de-paraffinized sections using eBioscience™ IHC Antigen Retrieval Solution - Low pH (10X) (Product # 00-4955-58) diluted to 1X solution in water in a microwave at 100 degree Celsius for 10 minutes. Following antigen retrieval, the sections were blocked with 2% normal goat serum in 1X PBS for 45 minutes at room temperature and then probed with or without CD45 Polyclonal Antibody (Product # PA5-96061) at 1: 100 dilution in 0.1% normal goat serum overnight at 4 degree Celsius in a humidified chamber. Detection was performed using Goat anti-Rabbit IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor Plus 488 (Product # A32731) at a dilution of 1:2000 in 0.1% normal goat serum for 45 minutes at room temperature. ReadyProbes™ Tissue Autofluorescence Quenching Kit (Product # R37630) was used to quench autofluorescence from the tissues. Nuclei were stained with DAPI (Product # D1306) and the sections were mounted using ProLong™ Glass Antifade Mountant (Product # P36984). The images were captured on EVOS™ M7000 Imaging System (Product # AMF7000) at 20X magnification and deconvoluted externally.

View more figures on thermofisher.cn

□1 Reference

Flow Cytometry (1)

iScience

A combination therapy of bortezomib, CXCR4 inhibitor, and checkpoint inhibitor is effective in cholangiocarcinoma *in vivo*.

"PA5-96061 was used in Flow cytometry/Cell sorting to find that the inhibition of CXCR4 allowed bortezomib to sensitize CCA to anti-PD1 treatment, with a significant reduction of tumor burden and long-term overall survival."

Authors: Li L,Zhou Y,Zhang Y,Hu H,Mao HQ,Selaru FM

Year 2023

Species Human

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPERSS OR IMPLED WARRANTIES, EXPERSS OR IMPLED WARRANTIES OF MERCHANTABILITY, ETTINESS FOR ANY PARTICULAR PURPOSE, OR NON HERIOGENETY.

BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.