Performance guarenteed

# Chicken anti-Goat IgG (H+L) Cross-Adsorbed Secondary Antibody, Alexa Fluor™ 594

#### **Product Details**

Size	1 mg
Species Reactivity	Goat
Host/Isotype	Chicken / IgY
Class	Polyclonal
Туре	Secondary Antibody
Conjugate	Alexa Fluor™ 594
Excitation/Emission Max	590/618 nm
Immunogen	Gamma Immunoglobins Heavy and Light chains
Form	liquid
Concentration	2 mg/mL
Purification	purified
Storage buffer	PBS, pH 7.5
Contains	5mM sodium azide
Storage conditions	4° C, store in dark
RRID	AB_2535871
KKID	AD_2333071

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	1-10 μg/mL	0 Publication
Immunohistochemistry (Paraffin) (IHC (P))	-	0 Publication
Immunohistochemistry (PFA fixed) (IHC (PFA))	-	0 Publication
Immunohistochemistry (Frozen) (IHC (F))	-	0 Publication
Immunocytochemistry (ICC/IF)	0.5-10 μg/mL	0 Publication
Miscellaneous PubMed (Misc)	-	0 Publication


### **Product Specific Information**

This Chicken anti-goat antibody reacts with IgG heavy chains and all classes of immunoglobulin light chains from goat. Chicken secondary antibodies have gained popularity because they demonstrate a lower level of nonspecific binding. Chicken antibodies lack a classic "Fc" domain and will not bind to protein A or protein G, nor will they bind to mammalian IgG Fc receptors.

Product will be shipped at Room Temperature.


Product Images For Chicken anti-Goat IgG (H+L) Cross-Adsorbed Secondary Antibody, Alexa Fluor™ 594

1

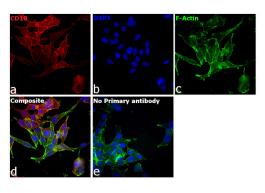


Goat IgG (H+L) Cross-Adsorbed Secondary Antibody (A-21468) in ICC/IF Immunofluorescence analysis of Chicken anti-Goat IgG (H+L) Cross-Adsorbed Secondary Antibody, Alexa Fluor™ 594 (Product # A-21468) was performed using LNCaP cells stained with CD10 Polyclonal Antibody (Product # PA5-47075). The cells were fixed with 4% paraformaldehyde for 10 minutes, permeabilized with 0.1% Triton<sup>™</sup> X-100 for 10 minutes, blocked with 2% BSA for 1 hour and labeled with 1:150 dilution of primary antibody for 3 hours at room temperature. The cells were probed with Chicken anti-Goat IgG (H+L) Cross-Adsorbed Secondary Antibody, Alexa Fluor™ 594 (Product # A-21468, 1:3000 dilution) in 0.1% BSA in PBS for 45 minutes at room temperature. Membrane localization of CD10 was seen in LNCaP (Panel a: Red). Nuclei (Panel b: Blue) were stained with Hoechst33342 (Product # H1399). F-actin was stained with Alexa Fluor® 488 Phalloidin (Product # A12379, 1:300) (Panel c: Green). Panel d represents the composite image. The specificity of the secondary antibody was proved by the absence of signal in HeLa (negative model for CD10) due to no primary antibody binding (Panel e). The images were captured at 40X magnification in CellInsight CX7 LZR High-Content Screening (HCS) Platform (Product # CX7A1110LZR) and externally deconvoluted (D.Sage et al./Methods 115 (2017) 28-41).

### Goat IgG (H+L) Cross-Adsorbed Secondary Antibody (A-21468) in ICC/IF

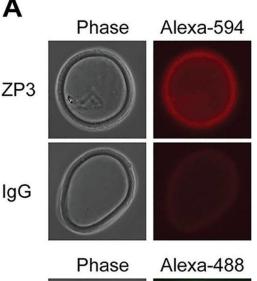

Amyloidogenic properties of the mouse ZP.A) The detection of amyloids in isolated mouse ZP was carried out using the amyloid conformation-dependent antibodies anti-fibrillar OC and anti-oligomer A11 in immunofluorescence analysis. Normal rabbit serum (NRS) served as a control. The anti-ZP3 antibody was used as a marker for the ZP with normal goat IgG serving as a control antibody. The corresponding phase images are shown for each fluorescent image. B) Intact ZP stained with 0.1% thioflavin S to detect amyloids. Scale bar =  $50 \ \mu\text{m}$ . C) ZP pellets stained with 0.2% Congo Red showed yellow-green birefringence (arrow) when examined under polarizing light and bright red fluorescence when examined with UV light. Scale bar =  $10 \ \mu\text{m}$ . Image collected and cropped by CiteAb from the following publication (https://pubmed.ncbi.nlm. nih.gov/26043223), licensed under a CC BY license.

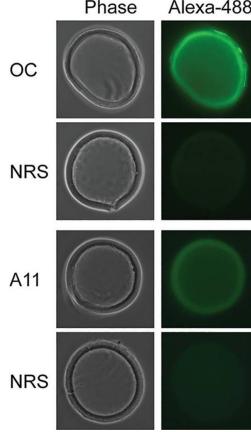

D	Ct	ΔPodA17	Ct diabetic	ΔPodA17 diabetic
ADAM17				
WT1				
MERGE			Er son Sala	and a second

**Goat IgG (H+L) Cross-Adsorbed Secondary Antibody (A-21468) in IHC (PFA)** Analysis of diabetic PodA17 mice model. (A) Blood glucose levels at different time points after STZ treatment. Twelve weeks after STZ administration, Ct and PodA17 non-diabetic and diabetic mice were analyzed for (B) 24-hour urinary albumin excretion, 24-hour urinary TNF- excretion, and serum TNF- level (n = 10 per group). (C) Glomerular ADAM17 activity in Ct and PodA17 non-diabetic and diabetic mice (n = 6 per group). (D) Immunofluorescence double staining showing expression of ADAM17 and WT1 in Ct and PodA17 mice, and their colocalization (merge). Note the increased staining for ADAM17 in the glomerulus of a Ct diabetic mice and its complete absence in podocytes of PodA17 mice, both diabetic and non-diabetic. Original magnification, X400 (\*p < 0.05, \*\*\*p 0.001; Student&quotes t test comparing diabetic mice, data are means  $\pm$  SEM) Image collected and cropped by CiteAb from the following publication (https://pubmed.ncbi.nlm.nih.gov/33634991), licensed under a CC BY license.

## **В 88 References**

Distinct effects of Fgf7 and Fgf10 on the terminal differentiation of murine bladder urothelium revealed using an organoid culture system. BMC Urol (2023)

Single-cell transcriptomic analysis of the identity and function of fibro/adipogenic progenitors in healthy and dystrophic muscle. iScience (2023)

SuPAR mediates viral response proteinuria by rapidly changing podocyte function. Nat Commun (2023)

Insulin-Like Growth Factor1 Preserves Gastric Pacemaker Cells and Motor Function in Aging via ERK1/2 Activation. Cell Mol Gastroenterol Hepatol (2023)

Cholesterol efflux pathways hinder KRAS-driven lung tumor progenitor cell expansion. Cell Stem Cell (2023)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and parts or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, this warranty as pecifications and dors accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty as one extend to anyone other than the Buyer. Any model or sample turnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. No OTHER WARRANTES, EXPRESS OR INNELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE EMEMPED FOR NON-CONFERNING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPARE, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION. THERE IS NO BELICATION TO REPAR. REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION. THERE IS NO BULGATION TO REPAR. REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS NURING THE WARRANT FOR OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER', (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product is intended for research only and is not to be used for any other purpose, including without limitation, unautherized commercial uses, in vitro diagnostic uses, or any type of consumption to human or